

syngo MammoReport

SP

Planning Guide

© Siemens AG 2005

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

English

Doc. Gen. Date: 10.05

Chapter	Page	Revision
All	All	01

Document revision level

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

Copyright notice

The DICOM functionalities in the *syngo* MammoReport system are based on the libraries from Merge Technologies Inc.

Celsius is a trademark of Fujitsu Siemens Computers.

Windows XP Professional is a trademark of Microsoft Corporation.

syngo is a trademarks of Siemens AG.

Coronis and MediCal Pro are trademarks of Barco.

Dome and CXtra are trademarks of Planar Systems, Inc.

Softcopy reading software by MeVis BreastCare,
copyright © 2001 - 2005 by MeVis BreastCare.



	Page
1 General Remarks	1-1
System configuration	1-1
Remarks	1-1
Safety	1-1
Requirements	1-1
System overview.	1-2
Network	1-4
Software installation	1-4
Software configuration.	1-4
Hardware configuration	1-4
Acronyms and abbreviations	1-5
2 Room Planning	2-1
Sample room plan and movement ranges.	2-1
Room planning example.	2-1
Permissible relative air humidity.	2-1
System: Front view	2-2
System: Side view.	2-3
3 Preparation for Installation	3-1
Uninterruptable Power Supply (UPS)	3-1
4 System Connections / Network	4-1
Overview of system connections	4-1
Cable length.	4-1
5 Technical Data	5-1
On-site electrical installation.	5-1
Heat dissipation data and environmental conditions	5-1
6 DICOM	6-1
DICOM compatibility.	6-1
7 Transport Conditions	7-1
Weight specifications	7-1
Component weights and dimensions	7-1
Table weight support	7-1
8 Network Nodes Checklist	8-1
syngo MammoReport Workstation	8-1
DICOM Conformance Statements.	8-1
DICOM Printer.	8-1

	Page
Acquisition workstation	8-2
CAD workstation	8-2
DICOM archive server / PACS	8-2
9 _____ Changes to previous version _____	9-1

System configuration

Remarks

- A note must be made on all room plans that are created by the Planning Departments which makes reference to the General Installation and Delivery Conditions of Siemens Medical Engineering. A copy of the General Installation and Delivery Conditions must accompany the plans when they are turned over to the customer.
- When planning the space, the documents which apply for the specific fire safety technical regulations must be observed.
- Unless otherwise stated, all dimensions are in mm.
- Room lighting:
For the constancy test according to DIN 6868-57, it must be possible to dim the lighting in rooms in which diagnostics is made on image display devices (monitors).
For the constancy test, it is also required that the value documented in the acceptance test must be exactly reproducible.
According to DIN 6868-57 (the international standard is still in preparation) the lighting in rooms where diagnosis is performed using image display devices (monitors) must meet the following requirements:
 - Adjustable, glare-free and reproducible setting of illuminance (e.g. dimmer with scale),
 - No reflections from windows, lamps and viewing boxes when the monitors are in their standard operating position.

Safety

The product-specific safety information contained in this document, as well as the general safety information must be observed.

NOTICE

Bear in mind that this product is intended to be used in a non-patient environment.

Requirements

It is very important that the customer is well prepared for the installation of the *syngo* MammoReport system. In order for the *syngo* MammoReport system to work as intended, an archive or printer solution is needed. Also make sure to complete the forms in this document as thoroughly as possible, and keep them available for the installation to facilitate the procedure.

NOTICE

VERY IMPORTANT NOTICE!

To enable monitor calibration according to the DIN6868-57 standard, SIEMENS strongly recommends that the SMfit ACT monitor calibration tools are purchased (including spot meter).

The DIN6868 plug-in which is also required is not included in SMfit ACT - it is available at a cost to the customer.

Instructions on how to perform quality control and calibrations using the SMfit ACT tool are described in the Quality Control Manual for *syngo* MammoReport.


System overview

The *syngo* MammoReport system's major parts are:

- 2 Displays (either CRT or TFT)
 - CRT displays: SIMOMED SMM21201P
 - TFT displays: DSB 2103-D-5MP, Barco MFGD 5621 HD , or Planar (model C5i)
- PC workstation
- Workflow keypad
- syngo compatible monitor
- syngo keyboard
- *syngo* MammoReport software
- Windows XP compatible printer, provided by the customer

NOTICE

The customer must provide a printer to be able to print reports.

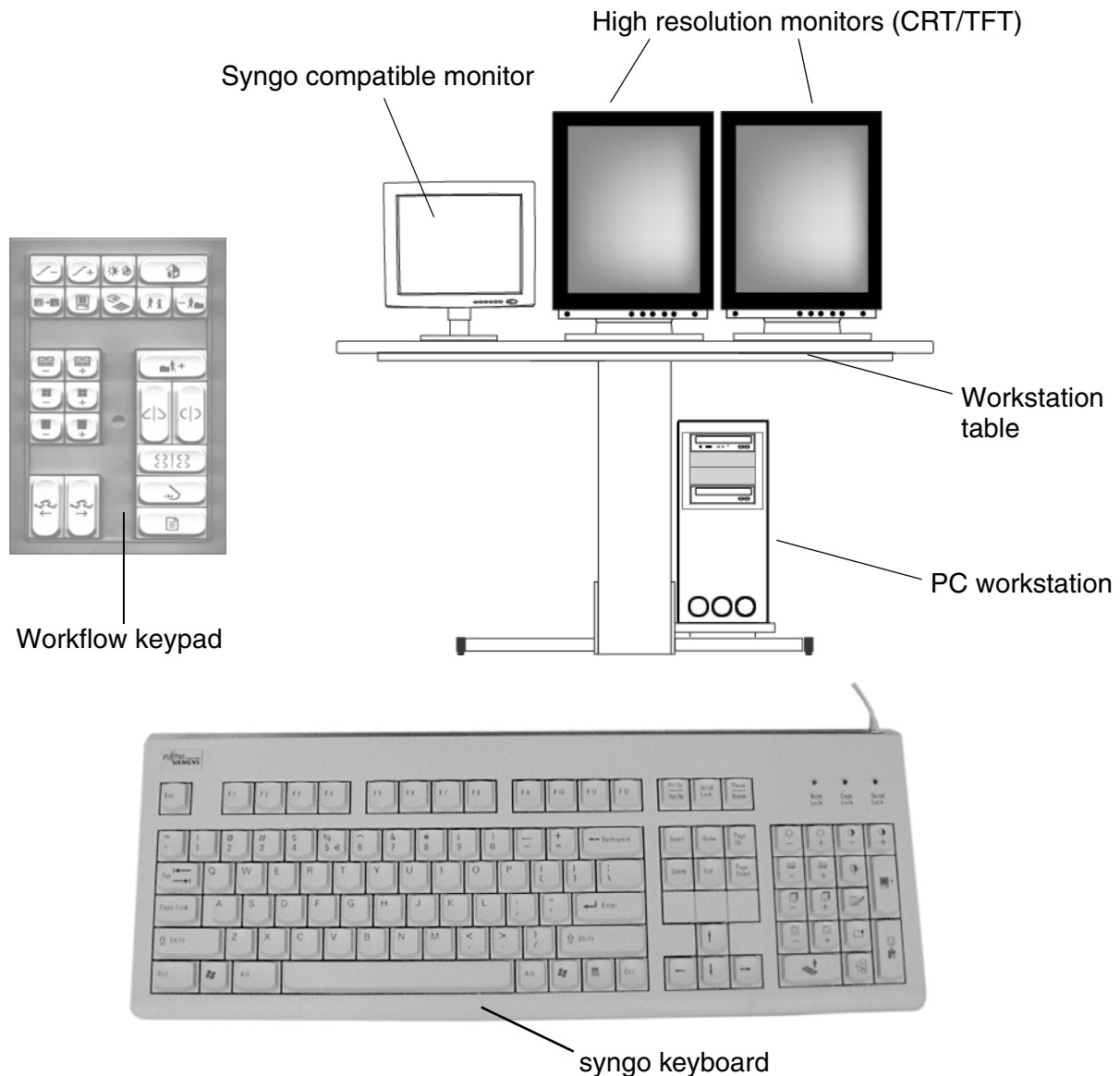

Additional options for the *syngo* MammoReport system are:

- 1 x 147 GB hard disk (for R610)
- 1 x 300 GB hard disk (for R630)

NOTE

For actual system options, see price book.

NOTE

The side table has to be provided by the customer.

If the workstation station table (option) is not ordered, the table must be provided by the customer.
--

Network

The *syngo* MammoReport workstation supports the TCP / IP protocol. Due to the data throughput that is required, we urgently recommend the use of a 100 Mb/s switched network or higher. Please note that any required network cabling, as well as the required connector outlets must be available on site. If necessary, an outside service contractor must be brought on site.

Software installation

The *syngo* MammoReport workstation is shipped with the softcopy reading software pre-installed, but will need customer specific settings configured.

Software configuration

The configuration of the *syngo* MammoReport workstation software requires approximately 2 hours.

Hardware configuration

The configuration of the *syngo* MammoReport workstation hardware, monitor calibration etc., requires approximately 2 hours.


Acronyms and abbreviations

AE Title (AET)	DICOM Application Entity Title. Must be unique in a DICOM network.
AWS	Acquisition Workstation.
Conformance	Each manufacture of a DICOM device has to provide such a statement which gives an overview of the product's DICOM capability.
CR	Computed Radiography IOD, e.g. used by the DIGISCAN M.
DICOM	Digital Imaging and Communication in Medicine.
DHCP	Dynamic Host Configuration Protocol. Reduces the complexity of configuring computers for TCP/IP networks. The following TCP/IP configuration can be dynamically assigned by a DHCP server: IP Address, Subnet masks, gateway and additional parameters like domain name.
DNS	The Domain Name Service is used to resolve the TCP/IP address from a "user friendly host name".
Hub	Distributor in a network using the star topology.
IP Address	32-bit address assigned to hosts using TCP/IP. An IP address belongs to one of four classes (A, B, C, D) and is written as 4 pairs separated with periods (decimal point format). The IP is part of the socket address. Here, IP means Internet Protocol.
SMPTE	Society of Motion Picture and Television Engineers
MG	Digital Mammography.
Q/R	Query / Retrieve is a service class which provides the ability to retrieve/transfer a well-identified set of images.
Port	Part of the socket address. Different DICOM services on the same host PC may use different port numbers.
Subnet mask	The subnet mask is used to extract network and sub network information from the IP address.
Router	A device that connects different networks. In a TCP/IP environment a router is used to connect computers that are located in different IP address ranges.
Service Class	A structured description of a service which is supported by cooperating DICOM Application Entities, e.g. Storage Service Class.
SCP	Service Class Provider, e.g. a Storage Service Class Provider provides the capability to receive DICOM images from a DICOM AET Storage SCU.
SCU	Service Class User, e.g. a DICOM Storage Service Class User can send an image to a DICOM Service Class Provider.
TCP	Transfer Control Protocol.
UI	User Interface.
USC	Uptime Service Center.
WINS	WINS provides a distributed database for registration and querying dynamic Net-BIOS names to IP addresses.

This page intentionally left blank.

Sample room plan and movement ranges

Room planning example

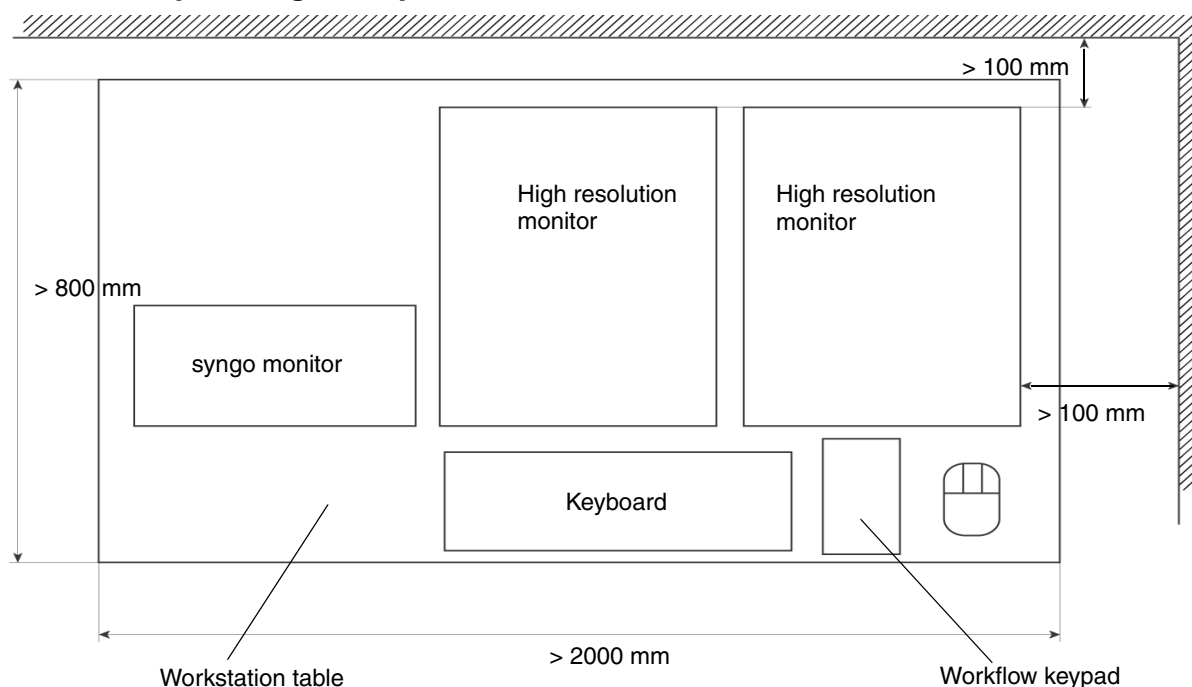

This example places the PC workstation underneath the workstation table. The room planning should also be prepared so that there are no power cables/transformers adjacent to the monitors. The monitors should also not be closer than 15 cm to other units.

To comply with the requirements in DIN 6868-57 and to provide optimal viewing conditions for diagnosis, lighting conditions in the room around the *syngo* MammoReport need to fulfil certain conditions. This is important to ensure that viewing of especially the darkest grayscale levels is not interfered with.

The *syngo* MammoReport monitors must be placed to avoid direct illumination and reflections on the monitors from light sources such as windows, lamps and light boxes. When viewing images on the *syngo* MammoReport, the illumination must come from diffuse light sources and it is recommended that the ambient light level is below 10 cd/m². The monitors should also be illuminated as evenly as possible from different directions.

NOTE

For cable lengths, see Page 4 - 1.


The *syngo* MammoReport should be installed in a room separate from the acquisition system.

The side table is provided by the customer.

Permissible relative air humidity


NOTE

Make sure that the room environment satisfies the requirements for relative air humidity, see "Heat dissipation data and environmental conditions" on Page 5 - 1. An air humidifier or an air dehumidifier might have to be installed.

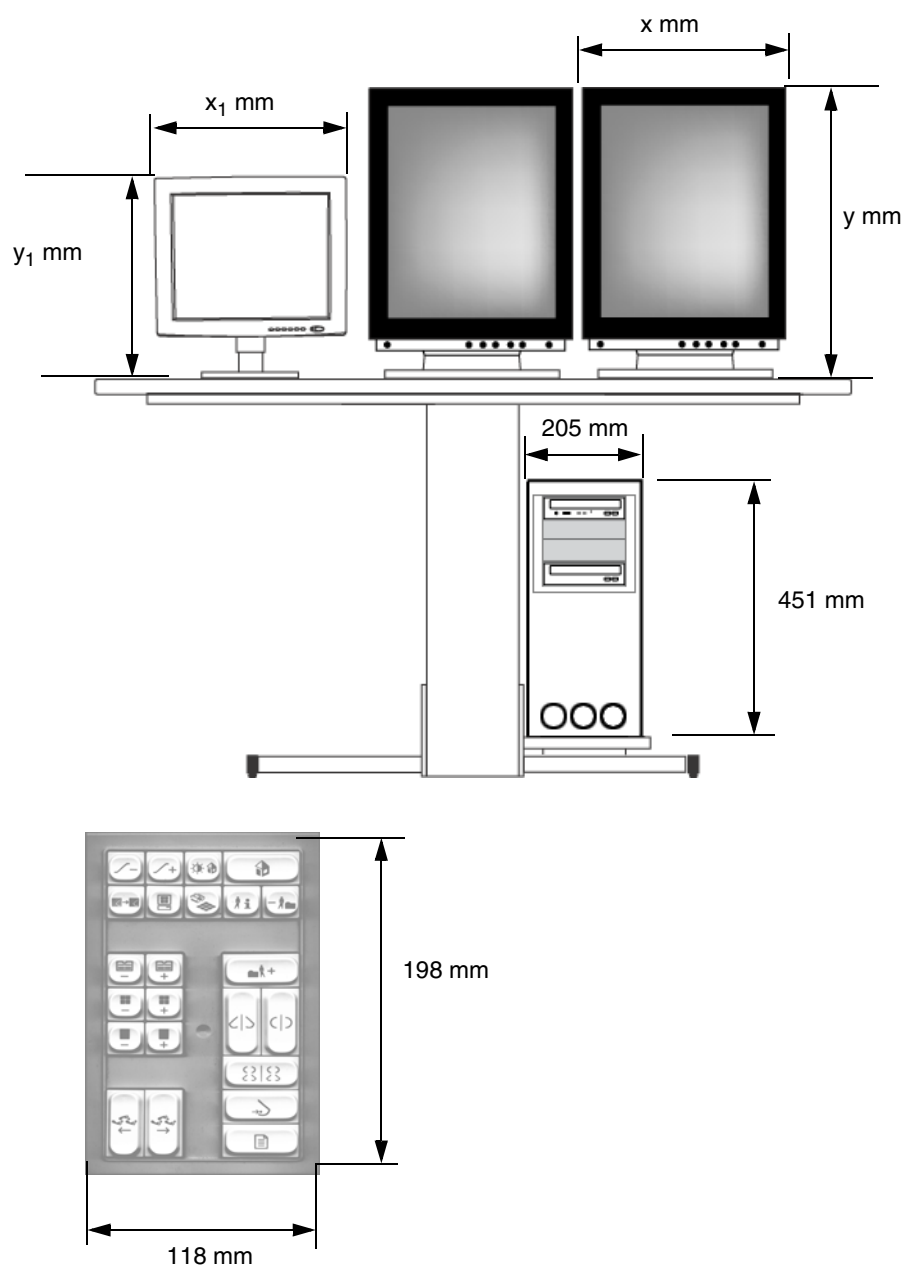System: Front view

	x resp. x_1	y resp. y_1
SIMOMED SMM21201P (CRT)	403 mm	580 mm
DSB 2103-D-5MP (TFT)	384 mm	521 mm
Planar C5i (TFT)	406 mm	512-590 mm (height adjustable)
BARCO MFGD 5621HD (TFT)	408 mm	592 mm
DSC 1904-D (syngo monitor)	420 mm	408-468 mm (height adjustable)

System: Side view

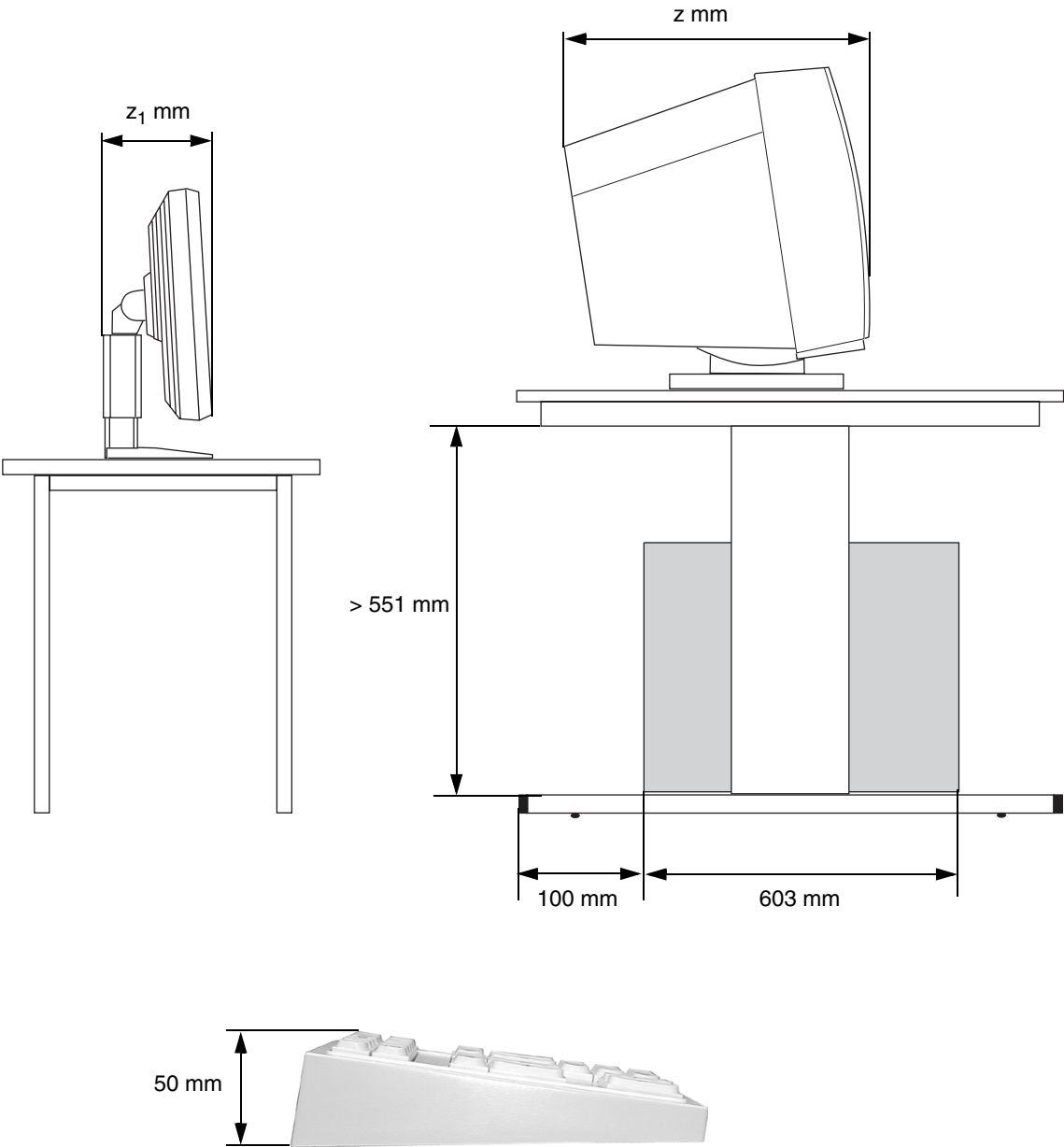

s

	z resp. z_1
SIMOMED SMM21201P (CRT)	520 mm
DSB 2103-D-5MP (TFT)	263 mm
Planar C5i (TFT)	239 mm
BARCO MFGD 5621HD (TFT)	217 mm
DSC 1904-D (syngo monitor)	210 mm

This page intentionally left blank.

NOTE

If additional SIENET components are added or if there is integration into a SIENET network, SHS/PACS P (Project Engineering) must be contacted.


Uninterruptable Power Supply (UPS)

CAUTION

To ensure that no data (patient data, image data) get lost if there is an unstable power supply and to protect the operating system (with all components), we recommend that you install a UPS.

This page intentionally left blank.

Overview of system connections

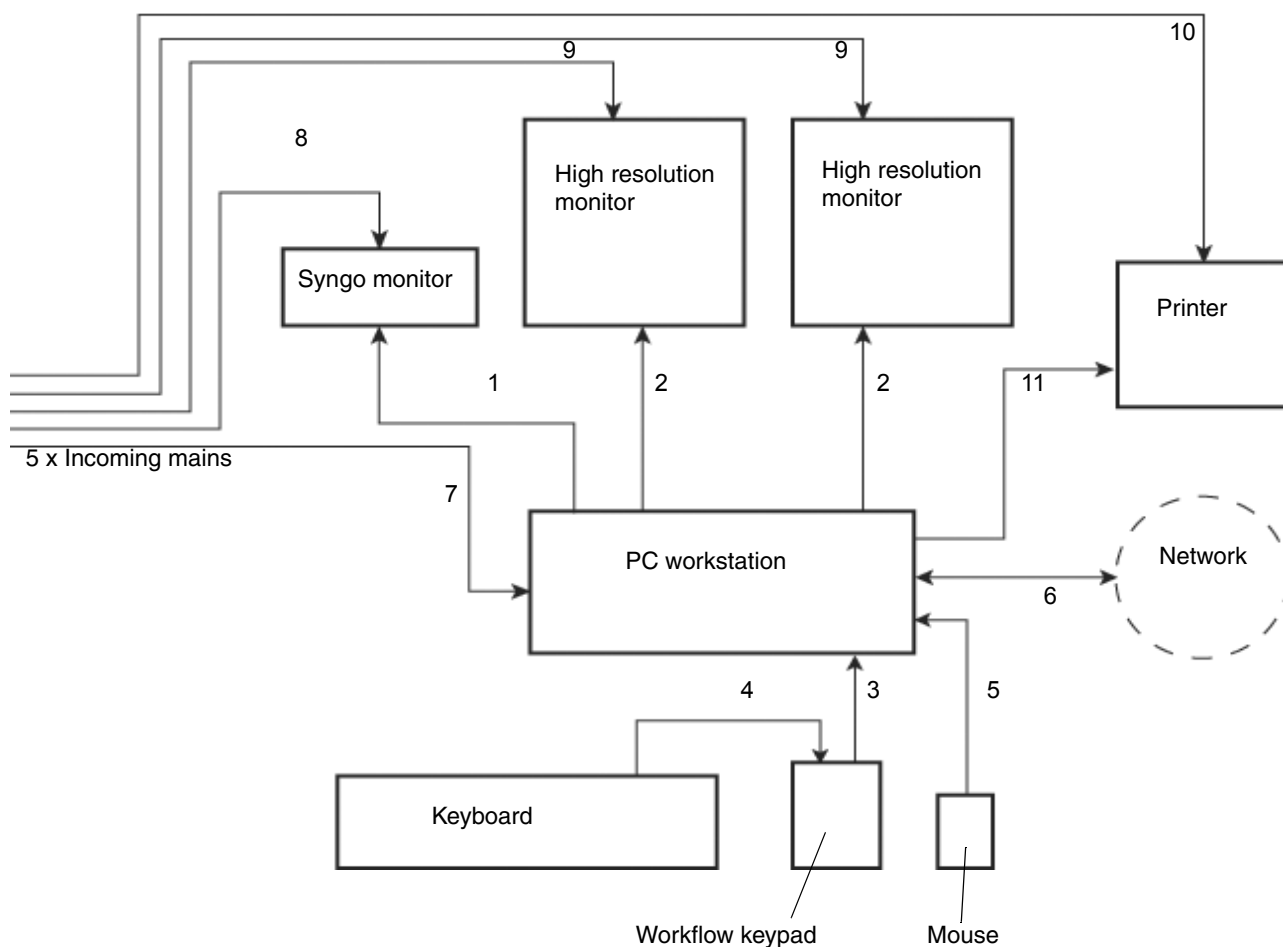

Cable length

	Cable between		Length
1	PC workstation	Syngo monitor	3 m
2	PC workstation	High resolution monitor	2 x 2 m
3	PC workstation	Workflow keypad	1 m
4	Workflow keypad	Keyboard	1 m
5	PC workstation	Mouse	1 m
6	PC workstation	Network	30 m
7	PC workstation	Incoming mains	2 m
8	Syngo monitor		3 m (DVI plus Power cable extension)
9	High resolution monitor		2 x 1.5 m
10	Printer		ref. to printer supplier
11	PC workstation	Printer	ref. to printer supplier

This page intentionally left blank.

On-site electrical installation

Proposal for building power distribution according DIN VDE 0107 or national regulations:

	Line voltage	Power line	Current consumption
PC workstation	100 - 240 V	AC 50/60 Hz	8.6 A @ 115 VAC, 4 A @ 230 VAC
syngo monitor (Siemens)	100 - 240 V	AC 50/60 Hz	1.0 A @ 100 VAC, 0.5 A @ 240 VAC
SIMOMED SMM21201P (CRT)	100 - 240 V	AC 50/60 Hz	2 A @ 100 VAC, 0.8 A @ 240 VAC
Planar C5i (TFT)	100 - 240 V	AC 50/60 Hz	2 A @ 100 VAC, 0.8 A @ 240 VAC
DSB 2103-D-5MP (TFT)	100 - 240 V	AC 50/60 Hz	2 A @ 100 VAC, 0.8 A @ 240 VAC
BARCO MFGD 5621HD (TFT)	90 to 264 V	AC 45/65Hz	2 A @ 100 VAC, 0.8 A @ 240 VAC

Heat dissipation data and environmental conditions

	Heat dissipation [kW]	Operation permissible ambient temperature	Operation permissible relative air humidity	Transport and storage permissible ambient temperature	Transport and storage permissible relative air humidity
PC workstation	<0.7	+10- +35°C	20-80 RH%	-30- +60°C	10-95 RH%
syngo monitor (Siemens)	<0.08	+ 5- +35°C	max. 75 RH% at 25 °C	-20- +60°C	max. 75 RH% at 25 °C
SIMOMED SMM21201P (CRT)	<0.2	+10- +35°C	20-80 RH%	-25- +70°C	20-80 RH%
Planar C5i (TFT)	<0.13	5- +35°C	max. 75 RH% at 25 °C	-20- +60°C	max. 75 RH% at 25 °C
DSB 2103-D-5MP (TFT)	<0.13	5- +35°C	max. 75 RH% at 25 °C	-20- +60°C	max. 75 RH% at 25 °C
BARCO MFGD 5621HD (TFT)	<0.13	0- +40°C	max. 75 RH% at 25 °C	-20- +60°C	max. 75 RH% at 25 °C

This page intentionally left blank.

DICOM compatibility

The following table gives a short overview of the products compatible with the acquisition workstation. For a more detailed information please check the DICOM conformance statements. The inter-operability with foreign products might be complex and must be clarified before starting the installation.

NOTE

For further information please refer to:

[http://www-td.med.siemens.de/Product Information/Planning/General...](http://www-td.med.siemens.de/Product%20Information/Planning/General...)

DICOM Planning Guide, TDIT-000.891.02... .

	Storage	Q/R	Print	Verification
<i>syngo</i> MammoReport Workstation	SCP	SCU	SCU	SCP/ SCU

This page intentionally left blank.

Weight specifications

The system will be delivered on 2 standard pallets, weighing approximately a total of 130 kg (with SIMOMED SMM21201P) resp. 90 kg (with DSB 2103-D-5MP).

Component weights and dimensions

	Weight [kg]	Dimension (HxWxD) [mm]
PC workstation	25	451x205x603
syngo monitor	8.7	408-468 x 420 x 210
SIMOMED SMM21201P	2 x 33	580x403x520
DSB 2103-D-5MP	2 x 13.3	521x384x263
Planar C5i (TFT)	2 x 6.9	512-590 x 406 x 239
BARCO MGFGD 5621 HD (TFT)	2 x 13.5	592x408x217
Workflow keypad	-	50x118x198

CAUTION

The High resolution monitors are heavy.

Always handle the high resolution monitors with two persons in order to avoid injuries.

Table weight support

- The Workstation table must be able to support a weight according to the Weight specifications.
- The side table must be able to support a weight of 6 kg resp. 10 kg minimum (the syngo monitor).

CAUTION

The supported weight specified for the tables is not valid if other equipment not belonging to the *syngo* MammoReport system is placed on the tables.

This page intentionally left blank.

syngo MammoReport Workstation

Parameter	Configuration	Remarks
Local Host (computer name)		Host name of the <i>syngo</i> MammoReport Workstation. Do not use more than 10 characters and no special characters
IP address Subnet mask		Host address and subnet mask
IP Gateway(s) address		IP addresses of the router
Own DICOM AE title:		The AET name must be unique in the network. Usually the AE title is set to the host name (upper case is recommended)
Port no.		
SCR AET SCR Port no.		Only used for CAD SR

Tab. 1 *syngo* MammoReport workstation data sheet.

DICOM Conformance Statements

The necessary DICOM conformance statements must be available. Check the Interoperability Database for additional information:

URL: www.td.med.siemens.de. Collection: Product Information/DICOM.

Did you check the Interoperability Database?

☐ Yes

☐ No

Did you check the Conformance Statements?

☐ Yes

☐ No

DICOM Printer

Parameter	Configuration	Remarks
Manufacturer		
Model		
Host name		
TCP/IP address		
Logical name		
DICOM AE title:		
Port no.		
Dmin		
Dmax		

Parameter	Configuration	Remarks
Film size [pixels]	Portrait: Landscape:	

Tab. 2 DICOM printers data sheet

Acquisition workstation

Parameter	Configuration	Remarks
Host name		
IP address		
DICOM AET		
Port no.		

Tab. 3 Acquisition workstation data sheet.

CAD workstation

Parameter	Configuration	Remarks
Host name		
IP address		
DICOM AET		
Port no.		

Tab. 4 CAD workstation data sheet.

DICOM archive server / PACS

Parameter	Configuration	Remarks
Location		
Archive manufacturer		
SW version		
Logical name		This is the name as it will appear in the DICOM Setup user interface
Host (Node) name		
IP address		Not necessary if a DNS or WINS name resolution is available

Tab. 5 DICOM server data sheet

Parameter	Configuration	Remarks
DICOM: General AET Port no.:		
Storage AET Port no.		If applicable
Query/Retrieve AET Port no.		If applicable
Query/Retrieve IM AET Port no.		If applicable
Verification Port no.		If applicable
Data model:		Patient root or study root

Tab. 5 DICOM server data sheet

This page intentionally left blank.

Chapter	Page	Change
all		initial version

This page intentionally left blank.